



EUD VALVE ACTUATOR/24 VDC

(AUTO SAMPLER-MALE CONNECTOR)

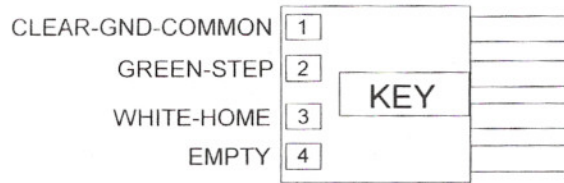


CABLE COMES WITH ACTUATOR

(GC-FEMALE CONNECTOR)

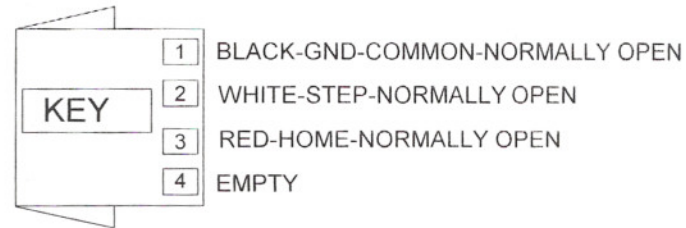


STANDARD 120V VALCO VALVE ACTUATOR

(AUTO SAMPLER-MALE CONNECTOR)

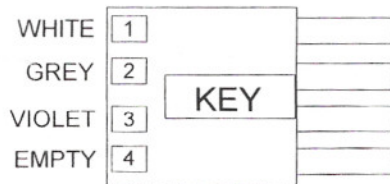

(GC-FEMALE CONNECTOR)

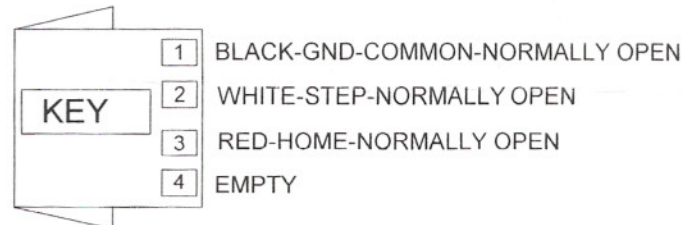

5030 AUTO SAMPLE CABLE

3' CABLE 3WIRE WITH CHASSIS GROUND


GC FEMALE CONNECTOR FOR RELAY

WHITE WIRE-RELAY#1 (N/O)-STEP

RED WIRE- RELAY#2 (N/O)-HOME

BLACK WIRE- COMMON JUMP RELAY 1&2


CHASSIS GROUND TO CHASSIS GROUND

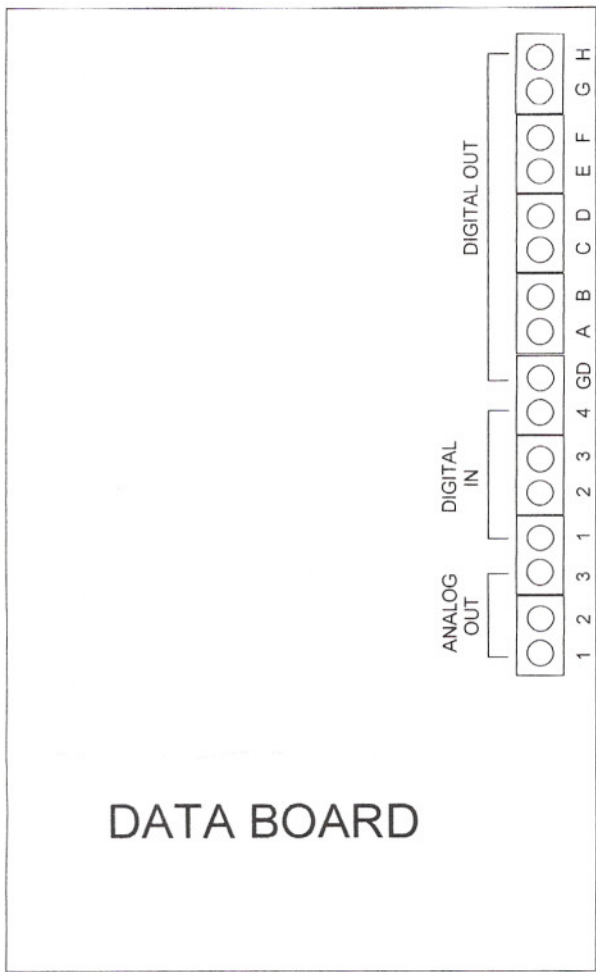

YEL/RED-(HOME) OR ANY OPEN EVENT

YEL/BRN-(STEP) WIRE OR ANY OPEN EVENT

RELAY BOARD

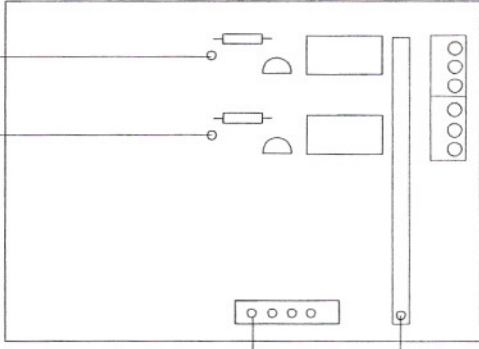

CHASSIS GND

12v BLU/RED WIRE TO AMP BOARD

